


Edwards ThruPort Systems EndoVent Pulmonary Catheter Kit

Product Overview

The EndoVent pulmonary catheter is indicated for use in patients undergoing cardiopulmonary bypass. It is intended to remove blood from the pulmonary artery, assist in decompressing the heart, maintain a dry operative field and monitor pulmonary artery pressure.

- Double lumen, 8.3 Fr catheter
- 5 cm interval markings to assist with placement
- It is pre-shaped and flexible to facilitate percutaneous placement through the internal jugular vein or the subclavian vein
- Integrated balloon to flow-direct the catheter into the pulmonary artery
- Protected by a contamination guard
- 3-way stopcock allows for pulmonary artery pressure monitoring

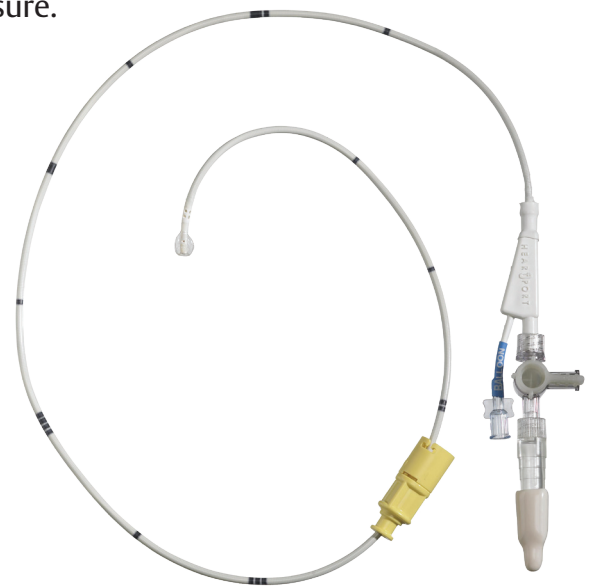

Model	Description	Size	Length
EV	EndoVent pulmonary catheter kit	8.3 Fr (2.8 mm)	25.6" (65 cm)


Contents:

- 1 EndoVent pulmonary catheter 8.3 Fr (2.8 mm)
- 1 Syringe (1.5 ml)
- 1 Catheter introducer sheath kit 9 Fr (3.0 mm)



Edwards

Edwards ThruPort Systems

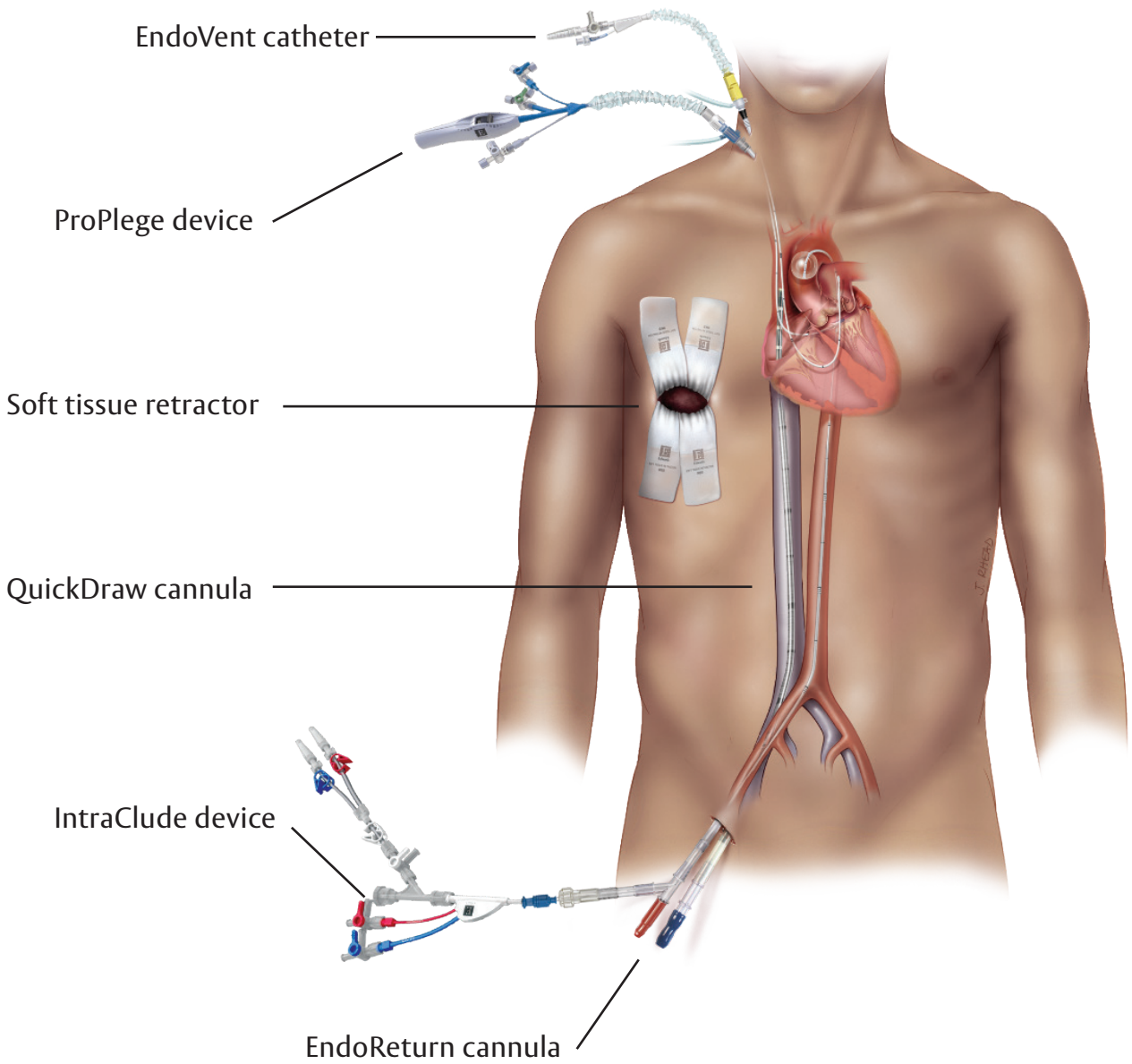

Caution: Federal (United States) law restricts this device to sale by or on the order of a physician. See instructions for use for full prescribing information, including indications, contraindications, warnings, precautions and adverse events.

Edwards, Edwards Lifesciences, the stylized E logo, EndoReturn, EndoVent, IntraClude, ProPlege, QuickDraw, and ThruPort are trademarks of Edwards Lifesciences Corporation. All other trademarks are the property of their respective owners.

© 2017 Edwards Lifesciences Corporation. All rights reserved. PP--US-1961v1.0

Edwards Lifesciences • One Edwards Way, Irvine CA 92614 USA • edwards.com



Edwards