

Temporary Transvenous

Pacing Catheters

Pacing, on demand

Whether you require temporary or instantaneous pacing, or hemodynamic monitoring with the security of pacing, Edwards Lifesciences will meet your needs. Each pacing catheter is part of the Swan-Ganz pulmonary artery catheter family and is backed by technical support, clinical and interactive training, educational tools and more.

Bipolar pacing catheter


- Designed for temporary right ventricular endocardial pacing when the patient's condition warrants
- Balloon tip facilitates insertion by flow direction
- A pair of electrodes at the tip provides capabilities for bipolar pacing



Swan-Ganz Paceport pulmonary artery catheter

The Paceport catheter facilitates temporary ventricular pacing and delivers basic hemodynamic monitoring, as well as:

- Right heart catheterization for right heart pressure measurements (PAS, PAD, PAOP) for diagnostic purposes
- Single point-in-time calculations of cardiac output using bolus thermodilution for diagnosing cardiac function
- Fluid infusion and RV pressure monitoring when the pacing probe is not inserted
- Available in oximetry Paceport catheter versions

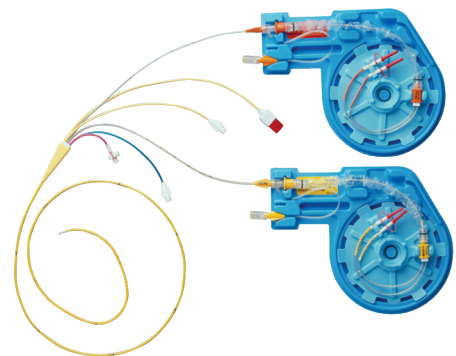

Chandler transluminal pacing for V-Pacing

- The Chandler transluminal V-Pacing probe (model D98100) can be used with multiple Paceport catheter models for standby ventricular pacing when the patient's condition warrants



Edwards

Pacing and thermodilution Paceport catheters

Model	Description	French Size	Length (cm)	Infusion Lumens
931F75	Paceport	7.5	110	3
D98100	Chandler transluminal V-Pacing probe	-	135	-
D200F7	Pacing thermodilution	7	110	1
D97120F5	Bipolar pacing catheters	5	90	1
D97130F5	Bipolar pacing (J-tip) catheters	5	90	1
780F75M	Oximetry Paceport	7.5	110	3

For more than 50 years, Edwards Lifesciences has been helping you make proactive clinical decisions and advance the care of surgical and acutely ill patients across the continuum of care.

To learn more visit www.Edwards.com/gb/SwanGanz

For professional use. For a listing of indications, contraindications, precautions, warnings, and potential adverse events, please refer to the Instructions for Use (consult eifu.edwards.com where applicable).

Edwards devices placed on the European market meeting the essential requirements referred to in Article 3 of the Medical Device Directive 93/42/EEC bear the CE marking of conformity.

Edwards, Edwards Lifesciences, the stylized E logo, Chandler, Paceport, Swan, and Swan-Ganz are trademarks of Edwards Lifesciences Corporation.

© 2021 Edwards Lifesciences Corporation. All rights reserved. PP--EU-2583 v2.0

Edwards Lifesciences • Route de L'Etraz 70, 1260 Nyon, Switzerland • edwards.com



Edwards