

Edwards General Surgery Catheters

A Spectrum of Solutions for Identification and Removal of Hepatic Calculi


Edwards Lifesciences offers a range of products for the identification and removal of hepatic calculi, including cholangiography and laparoscopic cholangiography catheters, Fogarty® biliary balloon probe, and Fogarty irrigation catheters.

Cholangiography Catheter

- Allows introduction of contrast media into cystic and common bile ducts
- Clear vinyl body makes it easy to detect air bubbles, which could be mistaken for stones
- Three-way luer lock hub enables the introduction of irrigation fluid and contrast media
- Black marking at the distal tip extends 1.5 cm proximally and facilitates placement
- Beveled tip facilitates introduction into common duct and reduces the risk of accidental withdrawal

Laparoscopic Cholangiography Catheter

- Cholangiography catheter also available in 80 cm length to meet the needs for laparoscopic procedures

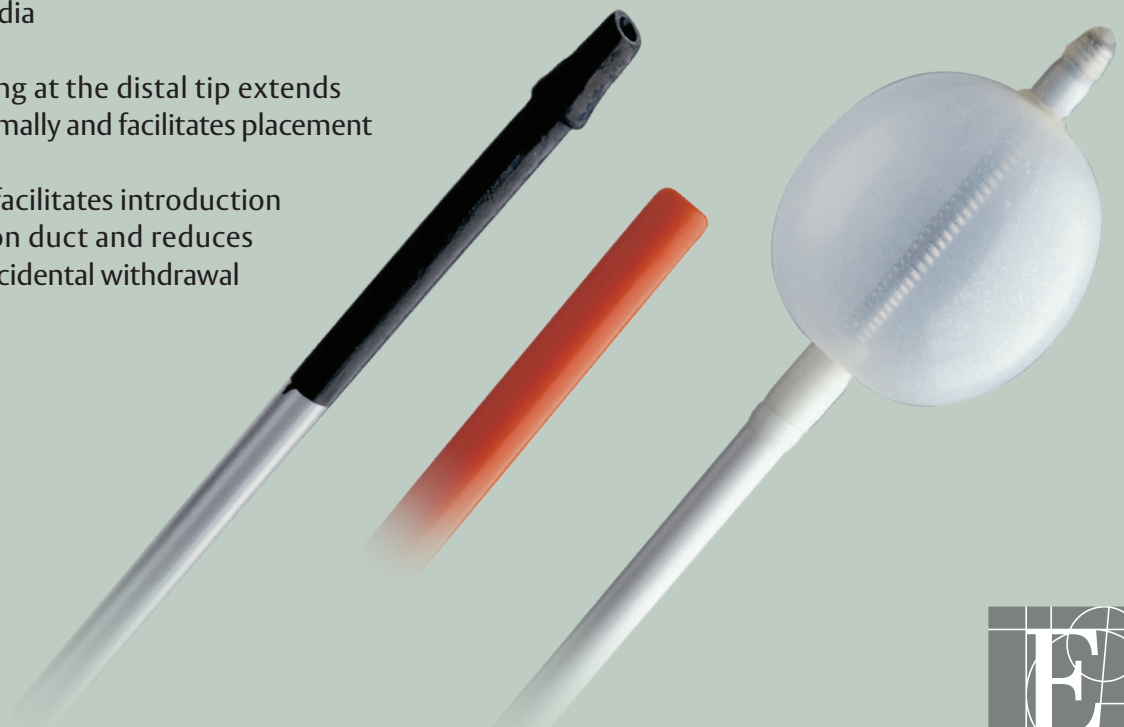

Edwards

Fogarty Biliary Balloon Probe

- Balloon catheter designed to atraumatically remove hepatic calculi
- Optimally thick balloon resists puncture by the calculi
- Radiopaque body is stiff enough to manipulate, yet flexible enough to minimize damage to the duct wall

Fogarty Irrigation Catheter

- 10F catheter designed with the Fogarty biliary balloon probe to irrigate the biliary tree after intraductal stone removal

SPECIFICATIONS			
Model Number	Catheter French Size	Length	Maximum Liquid Capacity
CHOLANGIOGRAPHY CATHETER			
420404F	4F	40 cm	
420406F	6F	40 cm	
LAPAROSCOPIC CHOLANGIOGRAPHY CATHETER			
420804F	4F	80 cm	
FOGARTY BILIARY BALLOON PROBES			
420235F	5F	23 cm	0.4 ml
410405F	5F	40 cm	0.4 ml
410236F	6F	23 cm	1.25 ml
FOGARTY IRRIGATION CATHETER			
220804F	4F	80 cm	
220806F	6F	80 cm	
2202310F	10F	23 cm	

For professional use. CAUTION: Federal (United States) law restricts this device to sale by or on the order of a physician. See instructions for use for full prescribing information, including indications, contraindications, warnings, precautions and adverse events.

Edwards Lifesciences devices placed on the European market meeting the essential requirements referred to in Article 3 of the Medical Device Directive 93/42/EEC bear the CE marking of conformity.

Edwards, Edwards Lifesciences, the stylized E logo, and Fogarty are trademarks of Edwards Lifesciences Corporation or its affiliates. All other trademarks are the property of their respective owners.

© 2016 Edwards Lifesciences Corporation. All rights reserved. PP- -US-0636 v1.0

Edwards Lifesciences • One Edwards Way, Irvine CA 92614 USA • edwards.com



Edwards